ASPIRATORE
ASPIRATOR
ASPIRATEUR
ABSAUGER
ASPIRADOR

ASKIR 36BR

30751/668 – Rev. 1 del 29.08.2012

Italian Medical Touch

CA-MI

MANUALE D'USO
INSTRUCTION MANUAL
MODE D'EMPLOI
HANDBUCH
MANUAL DE ISTRUCCIONES

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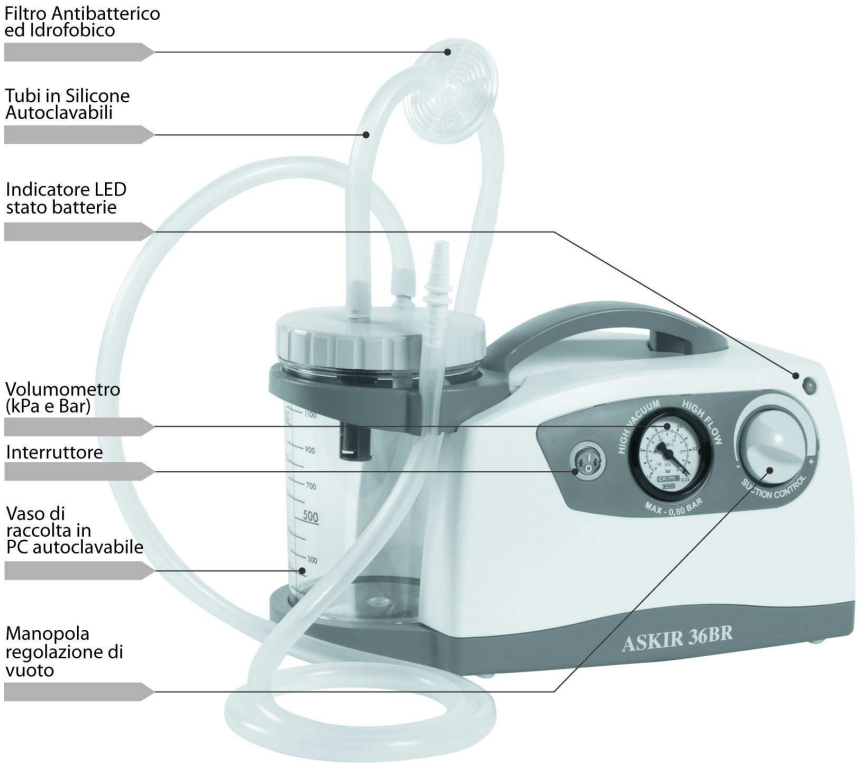
ASKIR 36BR is a surgical aspirator with the following electrica features : 14V ⎓ 4A with AC/DC adapter mod. UE60-140429SPA3 of FUHUA (input: 100-240V~ - 50/60Hz - 100VA) or Internally powered equipment (Pb Battery 12V ⎓ 4A) or with cigarette lighter adapter (12V ⎓ 4A). ASKIR 36BR is a desk-type electric suction unit for the aspiration of body liquids, oral, nasal and tracheal aspiration in adults or children. The appliance is designed for easy transport and almost continuous use.
The unit is operated by means of an electronic power management system, also monitoring the internal battery operation and status.
A LED system positioned on the front panel indicates appliance start-up and battery charge status.
The unit is equipped with suction regulator on the front panel and polycarbonate autoclavable jar complete with overflow valve.
Thanks to these characteristics and to its functions, this device is particularly suitable for different applications: utilization in hospital wards, for tracheotomy, for suction of body liquids and for minor surgery.

GENERAL WARNING

READ INSTRUCTION MANUAL CAREFULLY BEFORE USE
ONLY HIGHLY QUALIFIED STAFF USE RESERVED
THE INSTRUMENT MUST NOT BE DISASSEMBLED. FOR A TECHNICAL SERVICE ALWAYS CONTACT CA-MI

IMPORTANT SAFETY RULES

1. Check the condition of the unit before each use. The surface of the unit should be carefully inspected for visual damage. Check the mains cable and **don't connect to power** if damage is present.
2. Before connecting to power, always check the manufacture's label to ensure the electrical design of the unit and style of plug comply with the electrical network to be used.
3. If the plug supplied with the appliance is incompatible with the mains electricity socket, contact qualified staff for replacement of the plug with a suitable type. The use of simple or multiple and / or extension adapters is not generally recommended. Whenever their use is indispensable, use those in compliance with safety regulations, however paying attention not to exceed the maximum power supply limits, which are indicated on the adapters and extensions.
4. Pay special attention to the following:
 - Only use original accessories and components;
 - The device can be used only with the bacteriological filter;
 - Never immerge the appliance into water;
 - Position the appliance on flat stable surfaces;
 - Position the device in a way that the air inlets on the back aren't obstructed;
 - Never use the device in environments which have anaesthetic mixtures inflammable with air, oxygen or nitric oxide;
 - Don't touch the device with wet hands and always prevent the appliance coming into contact with liquids; Don't use the appliance if the plug or the AC/DC power supply are deteriorated or wet (send it promptly back to an authorised service center or the CA-MI technical service)
 - Keep off the reach of children or not capable people without supervision;
 - Don't leave the appliance connected to the power supply socket when not in use;
 - Hold plug to remove from socket. Don't pull on cable
 - Preserve and use the medical device in environments protected from atmospheric factors and at a distance from heat sources.
 - Don't use the device thoracic drainage.
5. For repairs, exclusively contact CA-MI technical service and request the use of original spare parts. Failure to comply with the above can jeopardise the safety of the device.
6. **This medical device must be destined exclusively for the use for witch it has been designed ad described in this manual.** Don't use for anything other than the use defined by the manufacturer. The manufacturer will not be responsible for damage due to improper use or connection to an electrical system not complying with current regulations.
7. The device is intended for use in the electromagnetic environment outlined in this manual. Install and use as directed.
8. Instrument and accessories discharging must be done following current law regulations in every country of use.
9. The lead battery integrated in the device is not to be considered as an ordinary domestic waste. Such a component must be disposed of in a specific collection centre in order to be recycled.
10. None of electric or mechanical parts have been designed to be repaired by customers or end-users. Don't open the device, do not mishandle the electric / mechanical parts. Always contact technical assistance
11. The use of the device in ambient condition other those indicated herein can seriously jeopardise the safety and the technical parameters of the appliance.

The manufactured cannot be held liable for accidental or indirect damages should the device be modified, repaired without authorization or should any of its component be damaged due to accident or misuse.

Any minimal modification / repair on the device voids the warranty and does not guarantee the compliance with the technical requirements provided by the MDD 93/42/EEC (and subsequent changes) and its normatives.

INSTRUCTIONS FOR THE CORRECT DISPOSAL OF THE PRODUCT FOR THE PURPOSE AND EFFECT OF EUROPEAN COMUNITY DIRECTIVE 2002/96/EC::

At the end of its useful life the product must not be disposed of with domestic waste.
It may be delivered to the relevant separate waste collection centres prepared by municipal administrations, or at retailers that supply this service. Disposing of the product separately prevents possible negative consequences for the environment and for health, deriving from inadequate disposal.
It also allows the recovery of materials of witch it's made up in order to obtain an important saving of energy and resources. In order to point out the obligation to dispose of the electromedical appliances separately, the product has the sign of the crossed mobile waste container applied to it.

TECHNICAL CHARACTERISTICS

Model	**ASKIR 36BR**
Typology (MDD 93/42/EEC)	Class IIa Medical device
UNI EN ISO 10079-1 Classification	HIGH VACUUM / HIGH FLOW
Power Feeding	14V ⎓ 4A with AC/DC adapter mod. UE60-140429SPA3 of FUHUA (input: 100-240V~ - 50/60Hz - 100VA) or Internally powered equipment (Pb Battery 12V ⎓ 4A) or with cigarette lighter adapter (12V ⎓ 4A)
Current Consumption	4.0A
Maximum Suction Pressure (without jar)	-80kPa (- 0.80 bar)
Minimum Suction Pressure (without jar)	Less -40kPa (-0.40 bar)
Maximum Suction Flow (without jar)	36 l /min
Weight	4.39 Kg
Insulation Class (when used with the AC/DC adapter mod. UE60-140429SPA3)	Class I
Insulation Class (when used with an Internal battery)	Internally Powered Equipment
Insulation Class (when used with a car cigarette lighter cable)	Class II
Size	357 x 193 x 215 mm
Battery Holding Time	60 minutes
Battery Time Charge	120 ÷ 150 minutes
Accurancy of Vacuum Indicator	± 5%
Working Condition	Room Temperature: 5 ÷ 35° C Room humidity percentage 30 ÷ 75% RH Altitude: 0 ÷ 2000m s.l.m
Conservation condition and Transport	Room Temperature: - 40 ÷ 70° C Room humidity Percentage: 10 ÷ 100% RH

SIMBOLOGY

CE 0123	CE marking in conformity with EC directive 93/42/EEC and subsequent changes
	General warnings and/or specifications
	Consult the instruction manual
	Keep in a cool, dry place
	Conservation temperature: -40 ÷ 70°C
	Type B equipment
~	Alternate Current
	Battery
	Direct Current
	Insulation Class II (only when connected to a car cigarette lighter cable)
Hz	Mains Frequency
I	ON
0	OFF
LOT	Lot Number
SN	Serial Number
REF	Identification device

ACCESSORIES SUPPLIED

DESCRIPTION	CODE
COMPLETE ASPIRATION JAR 1000cc	RE 210001
CONICAL FITTING	RE 210420
TUBES SET 8 mm x 14 mm	51100/01
HYDROPHOBIC AND ANTIBACTERIAL FILTER (+ n°1 spare parts)	SP 0121
AC/DC ADAPTER	SP 0208
EUROPEAN POWER SUPPLY CORD FOR AC/DC ADAPTER	SP 0020
CIGARETTE LIGTHER CABLE	SP 0007/02

The filter is produced with (PTFE) hydrophobic material to prevent fluids entering the pneumatic circuit. It should be changed immediately if it becomes wet or if there is any sign of contamination or discolouration. If should also be changed if the unit is used with a patient whose risk of contamination is unknown. **Don't use the suction unit without the protection filter**. If the suction unit is used in an emergency or in a patient where the risk of contamination is not know the filter must be changed after each use.
Available under request with different versions with complete jar 2000cc.

WARNING: The medical device is provided without a specific suction probe. If this device must be used with a specific suction probe, the end user is responsible for making sure it complies with the EN 10079-1 regulation.

WARNING: Suction tubes for insertion in the human body purchased separately from the machine should comply with ISO 10993-1 standards on material biocompatibility.

CLEANING OF ACCESSORIES

Washing and / or cleaning the autoclavable jar as to be carried out as follows:

- Wear protection gloves and apron (glasses and face mask if necessary) to avoid contact with contaminating substances;
- Disconnect the jar from the device
- Disconnect all tubes from the jar and the protection filter
- Empty and dispose of the content and of the suction catheter according to the laws in force in your country;
- Separate all parts of the cover (overflow valve, o-ring);

After disposing of disposable parts and disassembling the jar wash in running cold water and rinse thouroughly.
Then soak in warm water (temperature shall not exceed 60°C). Wash thouroughly and if necessary use a non-abrasive brush to remove incrustations. Rinse in running warm water and dry all parts with a soft cloth (non-abrasive).
The jar and the cover can be autoclaved by placing the parts into the autoclave and running one sterilization stem cycle at 121°C (1 bar relative pressure – 15min) making sure that the jar is positioned upsidedown.
Mechanical resistance of the jar is guaranteed up to 30 cycles of sterilization and cleaning at the indicated conditions (EN ISO 10079-1). Beyond this limit the physical-mechanical characteristics of the plastic may decrease and replacement of the part is therefore recommended.
After sterilization and cooling at environment temperature of the parts make sure that these are not damaged.
Assemble the jar as follows:

- Place the overflow valve into its seat in the cover (under VACUUM connector)
- Insert floating valve keeping the o-ring towards the opening of the cage
- Place the o-ring into its seat around the cover
- After completing assembling operations always make sue that cover seals perfectly to avoid vacuum leackages or liquid exit

The aspiration tubes can be sterilized on autoclave using a sterilization cycle at 120°C (1 bar relative pressure – 15min. The conical connector can be sterilized on autoclave using a sterilization cycle at 121°C (1 bar relative pressure – 15min.
The device is ready for a new employment now.

DO NOT WASH, STERILIZE OR PUT IN AUTOCLAVE THE ANTIBACTERIAL FILTER

CLEANING DEVICE

Use a soft dry cloth with not-abrasive and not – solvent detergents.

BEFORE STRATING CLEANING OPERATIONS, MAKE SURE THAT INTERNAL PARTS OF DEVICE DON'T COME INTO CONTACT WITH LIQUIDS AND THE AC/DC SWITCHING ADAPTER ARE DISCONNECTED.

DISPOSAL OF WASTE BATTERIES - (Directive 2006/66/EC)

This symbol on the battery or on the packaging indicates that the battery provided with this product shall not be treated as household waste. By ensuring these batteries are disposed of correctly, you will help prevent potentially negative consequences for the environment and human health which could otherwise be caused by inappropriate waste handling of the battery. The recycling of the materials will help to conserve natural resources.
At the end of their life hand the batteries over to the applicable collection points for the recycling of waste batteries.
For more detailed information about recycling of this product or battery, please contact your local Civic Office, your household waste disposal service or the shop where you purchased the product.

RULES FOR RETURNING AND REPAIRING

COMPLYING WITH THE NEW EUROPEAN RULES, CA-MI INDICATES THE IMPORTANT POINTS TO PROTECT INSTRUMENT AND OPERATORS HYGIENE. THESE RULES MUST BE RESPECTED IN ORDER TO GUARANTEE HYGIENE AND SAFETY TO ALL THE PEOPLE OPERATING WITH THE INSTRUMENT TO OBTAIN QUALITY AND WELL BEING.

Before returning an instrument for repair, the external surfaces and all accessoris **MUST** be carefully disinfected with a cloth soaked in methylated spirits or hypochlorite-based solution. If CA-MI finds instrument not suitable for repairing due to clear signs of internal or external contamination, the same will be returned to customer with specification of NOT REPAIRED INSTRUMENT, accompanied by an explanation letter. CA-MI will decide if contamination is due to bad functioning or misuse. If contamination is due to bad functioning, CA-MI will substitute the instrument, only if a SALE RECEIPT and STAMPED GUARANTEE accompany the same.
CA-MI is not responsable for contaminated accessories, they will be substitute at customer's expenses.
For this reson it is **COMPULSORY** to carefully disinfect the external part of the instrument and accessories with a cloth soaked in methylated spirits or hypochlorite-based solutions. Put the instrument and accessories in a bag with indication of disinfecting.
We also request to specify the kind of fault, in order to speed up repairing procedures.
To this end, please read the instructions carefully in order to avoid damaging the equipment through improper use.
Always specify the fault encountered so that CA-MI can establish whether it falls into the category of the faults covered by the guarantee.

Guidance and manufacturer's declaration – Electromagnetic Emissions		
The surgical aspirator **ASKIR 36BR** is intended for use in the electromagnetic environment specified below. The customers or the user of the surgical aspirator **ASKIR 36BR** should assure that it's used in such an environment.		
Emissions Test	**Compliance**	**Electromagnetic environment - guidance**
Irradiated / Conducted emissions CISPR11	Group 1	The surgical aspirator ASKIR 36BR only used RF energy only for its internal functioning. Therefore its RF emissions are very low and are not cause interference in proximity of any Electronic appliances.
Irradiated / Conducted emissions CISPR11	Class [B]	The surgical aspirator ASKIR 36BR can be used in all environments, including domestic and those connected directly to the public mains distribution that supplies power to environments used for domestic scopes or environments feeds to you from batteries.
Harmonic emissions IEC/EN 61000-3-2	Class [A]	
Voltage fluctuations / flicker emissions IEC/EN 61000-3-3	Complies	

Guidance and manufacturer's declaration – Immunity Emissions			
The surgical aspirator **ASKIR 36BR** is intended for use in the electromagnetic environment specified below. The customers or the user of the surgical aspirator **ASKIR 36BR** should assure that it's used in such an environment.			
Immunity Test	**Level indicated by the IEC 60601-1-2**	**Compliance Level**	**Electromagnetic environments - guidance**
Electrostatic discharge (ESD) IEC/EN 61000-4-2	± 6kV on contact ± 8kV in air	The device doesn't change its state	Floors should be wood, conceret or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient / burst IEC/EN 61000-4-4	± 2kV power supply lines ± 1kV for input / output lines	The device doesn't change its state	Mains power quality should be that of a typical commercial environment or hospital or environments feeds to you from batteries.
Surge IEC/EN 61000-4-5	± 1kV differential mode	The device doesn't change its state	Mains power quality should be that of a typical commercial environment or hospital or environments feeds to you from batteries.
Loss of voltage, brief voltage interruptions and variations IEC/EN 61000-4-11	5%U_T (>95% dip U_T) for 0.5 cycle 40%U_T (>60% dip U_T) for 5 cycle 70%U_T (>30% dip U_T) for 25 cycle <5%U_T (>95% dip U_T) for 5 sec	-	Mains power quality should be that of a typical commercial environment or hospital If the user of the surgical aspirator ASKIR 36BR request that the appliance operates continuosly, the use of a continuity unit is recommended.
Magnetic field IEC/EN 61000-4-8	3A/m	The device doesn't change its state	The power frequency magnetic field should be measured in the intended installation location to assure that it's sufficiently low.
Note U_T is the value of the power supply voltage			

This section contains information regarding the conformity of the compliance with the IEC 60601-1-2 Standard.
The ASKIR 36BR surgical aspirator is an electro-medical device that requires particular precautions regarding electro-magnetic compatibility and which must be installed and commissioned according to the electro-magnetic compatibility information supplied.
Mobile and portable RF communication appliances (mobile phones, transceivers, etc..) can affect the medical system.
The use of accessories, transducers and cables different to those specified, with the exception of transducers and cables sold by the appliance and system manufacturer as spare parts, can lead to an increase in emissions or in a decrease of the immunity of the device or system.

Guidance and manufacturer's declaration – Immunity Emissions			
The surgical aspirator **ASKIR 36BR** is intended for use in the electromagnetic environment specified below. The customers or the user of the surgical aspirator **ASKIR 36BR** should assure that it's used in such an environment.			
Immunity Test	**Level indicated by the IEC 60601-1-2**	**Compliance level**	**Electromagnetic environments - guidance**
Conducted Immunity IEC / EN 61000-4-6	3Vrms 150kHz to 80Mhz (for non life-supporting devices)	$V_1 = 3$ V rms	The portable and mobile RF communication devices, including cables, must not be used closer to the ASKIR 36BR device, than the separation distance calculated by the equation applicable to the transmitter frequency. Recommended separation distance $d = [3.5 / V_1] \sqrt{P}$ $d = [12 / E_1] \sqrt{P}$ from 80 MHz to 800MHz $d = [23 / E_1] \sqrt{P}$ from 800 MHz to 2.5 GHz Where P is the maximum nominal output voltage of the transmitter in Watt (W) depending on the manufacturer of the transmitter and the recommended separation distance in metres (m). The intensity of the field from the fixed RF transmitters, as determined by an electro-magnetic study of the site[a], could be lower than the level of conformity of each frequency interval [b]. It is possible to check for interference in proximity to devices identified by the following symbol: 
Radiated Immunity IEC / EN 61000-4-3	3V/m 80MHz to 2.5GHz (for non life-supporting devices)	$E_1 = 3$ V / m	

Note 1: At 80 MHz and 800 MHz the interval with the highest frequency is applied
Note 2: These guide lines may not be applicable in all situations. The electro-magnetic propagation is influenced by the absorption and by reflection from buildings, objects and people.
a) The field intensity for fixed transmitters such as the base stations for radiotelephones (mobile and cordless) and terrestrial mobile radio, amateur radio devices, radio AM and FM transmitters and TV transmitters can not be theoretically and accurately foreseen.
To establish an electro-magnetic environment generated by fixed RF transmitters, an electro-magnetic study of the site should be considered. If the field intensity measured in the place where the device will be used surpasses the above mentioned applicable level of conformity, the normal functioning of the device should be monitored. If abnormal performance arises, additional measures such as changing the device's direction or positioning may be necessary.
b) The field intensity on an interval frequency of 150 kHz to 80 MHz should be less than 3 V/m.

Recommended separation distance between portable and mobile radio-communication devices and the monitor

The ASKIR 36BR surgical aspirator is intended to operate in an electro-magnetic environment where RF irradiated interferences are under control. The client or operator of the ASKIR 36BR device can help prevent electro-magnetic interference by keeping a minimum distance between the portable and mobile RF communication devices (transmitters) and the ASKIR 36BR device, as recommended below, in relation to the radio-communication maximum output power.

Maximum nominal output power of the Transmitter W	**Separation distance from the frequency transmitter (m)**		
	150 kHz to 80 MHz $d = [3.5 / V_1] \sqrt{P}$	80 MHz to 800 MHz $d = [12/E_1] \sqrt{P}$	800 MHz to 2.5 GHz $d = [23/E_1] \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters with a maximum nominal output power not shown above, the recommended separation distance in metres (m) can be calculated using the equation applicable to the transmitter frequency, where P is the maximum nominal output power of the transmitter in Watt (W) depending on the transmitter's manufacturer.

Note 1: At 80 MHz and 800 MHz the interval with the highest frequency is applied

Note 2: These guide lines may not be applicable in all situations. The electro-magnetic propagation is influenced by the absorption and by the reflection from buildings, objects and people.

MAINTENANCE

The **ASKIR 36BR** suction equipment does not need maintenance or lubrication.
It is, however to inspect the unit before each use. Unpack the instrument and **always check** integrity of plastic parts and AC/DC switching adapter , feeding cable, they might have been damaged during previous use. Connect cable to electrical network and turn switch on. Close the aspirator outlet with your finger and with suction regulator in maximum vacuum position check that the vacuum indicators reaches – 80kPa (-0.80 bar) minimum (internal battery). Rotate the knob from right to left and check the aspiration regulating control. The vacuum indicator should go down – 25kPa (-0.25 bar).
Verify that loud noises are not present, these can indicate wrong functioning. The device is protected by a safety fuse (**F 10A L 250V**) situated in the cigarette lighter cable. When replacing, always check the type and value as indicated.
Internally, the device is protected (see electrical specifications) by two fuses (**T 15A L 125V**) that cannot be reached from the outside. Therefore, contact the manufacturer to request the assistance of an authorized and qualified technician when they need to be replaced. If it's replaced make sure that its replacement is always the same type and value, as indicated.
The device is made up of a lead battery which cannot be accessed by outside. In order to replace it, consult the technical staff authorised by the manufacturer.

Fault type	Cause	Solution
1. Red light on	Battery run down	Hook up the power cord to the electricity mains, positioning the equipment power switch on 0.
2. No light	Defective AC/DC adapter or technical internal problem	Contact the technical service
3. No aspiration	Jar Cap badly screwed down	Unscrewed the cap, then rescrew it correctly
4. No aspiration	Lid seal not in its seat	Unscrew the cap and insert the seal properly in its seat
5. The float doesn't close	If the cap has been washed, ensure that the float is not partially detached	Insert the float into it's place
6. The float doesn't close	The float it's covered by dirty material	Unscrewed the cap, leave the and put in on autoclave
7. Low suction	Foam inside the jar	Fill the jar to 1/3 full of ordinary water
8. No aspiration due to flow leakage of mucus	Filter blocked	Replace filter
9. The Vacuum power on the patient side is either very low or absent	• Vacuum regulator set to minimum • Protection filter blocked or damaged • Connection tubes blocked, kinked or disconnected • Shut-off valve blocked or damaged • Pump motor damaged	• Turn the vacuum regulator clockwise and check the value of the vacuum on the gauge • Replace the filter • Replace or reconnect the tubes, check the jar connections • Empty the jar, or disconnect the tube from the jar and unblock the shut-off valve. The unit twill only work in the upright position • Contact the technical service
10. Noisy	Technical internal problem	Contcat the technical service
Faults 1 - 2 - 3 - 4 - 5 - 6 - 7 - 8 – 9 - 10	**None of the remedies has achieved the desired results**	**Contact the seller or CA-MI After-sales Assistance Service**

If the overfill security system it's activated, don't proceede with the liquid aspiration.
If the overfill security system doesn't work there are two cases:
1° case – If the overfill security system doesn't work the aspiration will be stopped by the bacteriological filter who avoid the liquid penetration inside the device.
2° case – If both the security system doesn't work, there is the possibility that liquid comes inside the device, in this case return the device to CA-MI technical service.

CA-MI Srl will provide upon request electric diagrams, components list, descriptions, setting instructions and any other information that can help the technical assistance staff for product repair.

CA-MI Srl cannot be held liable for accidental or indirect damages should the device be modified, repaired without authorization or should any of its component be damaged due to accident or misuse.

Any minimal modification / repair on the device voids the warranty and does not guarantee the compliance with the technical requirements provided by the MDD 93/42/EEC Directive and its normatives.


INSTRUCTIONS

Operation using with AC/DC switching adapter

- Place the unit on a flat, horizontal surface.
- Connect the short silicon tube, with antibacterial filter, to the suction connector. The other tube, with one end connected to the antibacterial filter with the other end to jar's lid connector where has been fixed the red float. When the 90% of the volume of the jar is reached there is the activation of the security float (the float close the aspiration connector on the jar) to the avoid liquid penetration inside the device.
- Connect the long silicon tube to the other jar's lid connector
- Connect the other end of the long silicon tube to the probe plastic connector then connect the suction probe to it.
- Connect the switching adapter to the device with the appropriate connector and insert the power cable plug to the power socket. To start the treatment press the **I** switch to turn it on
- Set the desired vaccum value (Bar / kPa) with the appropriate vacuum regulator. Turning the handle clockwise increase the vacuum value: these values can be read on the "vacuum indicator" instrument.
- To suspend and / o r terminate the treatment, press the switch againg and pull the plug out from the power socket.
- Unscrew the jar's lid and fill the jar 1/3 full or ordinary water (this for an easy cleaning operations and an rapid reaching of the functionally vacuum) then rescrew the lid on the jar correctly.
- To extract the accessories and start with cleaning.

WARNING: Ensure that the IN marker on the filter is on the side facing the collection jar lid and fitted into the "VACUUM".
A wrong connection causes immediate destruction in case of contact with sucked liquids.

Filter assembling

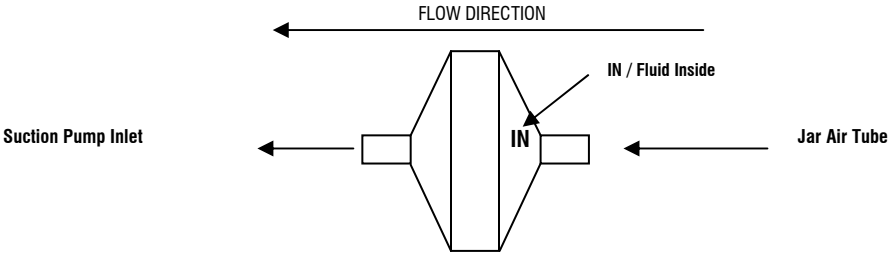If the state of the protective conductor of the installation, to which the medical device in question is connected, appear to be dubious, the use of equipment with a internal battery is recommended.

Operation using cigarette lighter DC 12V

- Connect the device's external plug 12V to the ligther plug with the cigarette lighter cable. Check the battery power status of the vehicle before the cigarette lighter cable. Press the switch to start suction
- Press the switch to the **I** position to turn it on.

WARNING: **Only use the originally supplied or recommended replacement cigarette ligther cables (view the chapter "Important Safety Rules")**

WARNING: Before using the device, check the battery power status. Before each use proceed with charging the battery. To maintain the device in good conditions, recharge the battery every 3 months (when not in use).

Recharging operations: to be able to charge the internal battery it is necessary to connect the universal switching adapter (mod. UE 60140429SPA3 of FUHUA) to the electric network for approx. 120 to 150 minutes with the main switch to position 0.
The battery's autonomy when fully charged is approx. 60 minutes with continued operations.


TAB. I – INDICATOR LIGHTS DURING OPERATIONS

When an external power supply (regardless of the state of the battery charger) and when the device is working (after having turned it on), the LED stays in a FIXED GREEN position.

LED Signal	Phase	Problem / Cause	Solution
Flashing Green Led	During rechanrge	Battery recharge running	Wait
Steady Greed Led	During recharge	Recharging cycle complete	Remove power supply
Steady Red Led	During battery operation	Flat battery	Start recharging cycle
Flashing Red Led	Device automatically turns off when the battery is flat	Battery completely flat	When the device is restarted the LED will flash red: begin the battery recharge cycle immediately
Steady Orange Led	During battery operation	Intermediate status	Guaranteed battery function / Recharge when the red LED signal comes on.

NEVER USE THE DEVICE WITHOUT JAR AND / OR PROTECTION FILTER

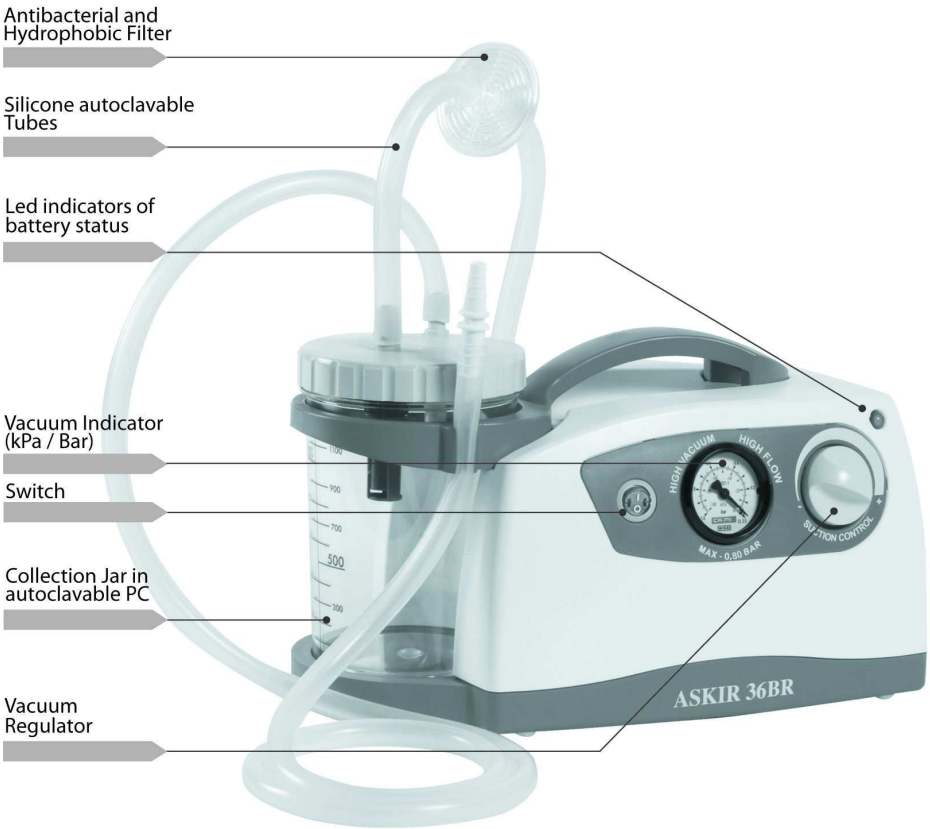MAKE SURE THAT CHILDREN AND/OR MENTALLY ILL PEOPLE DO NOT USE THE DEVICE WITHOUT ADULT SURVEILLANCE

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